


Autosampler Tubes and Caps

12 mL and 15 mL

AnalytiChem autosampler tubes are manufactured from virgin polypropylene (PP) with ultra-low metal content, making them ideal for applications involving ICP or ICP-MS analysis. These round bottom tubes are available in two convenient sizes: 16 x 100 mm (12 mL) and 17 x 100 mm (15 mL) and designed to meet the requirements of position racks used by most instrument manufacturers.

Every lot of autosampler tubes receives a Certificate of Volume analysis confirming the product's compliance with Class A Specification according to ASTM E1272 and ISO 384 per our ISO 17025 accredited procedure.

In addition, autosampler tubes are now available with a Trace Metals Certificate, providing verified results for over 65 elements to ensure ultra-low contamination and compliance for critical trace metal analysis.

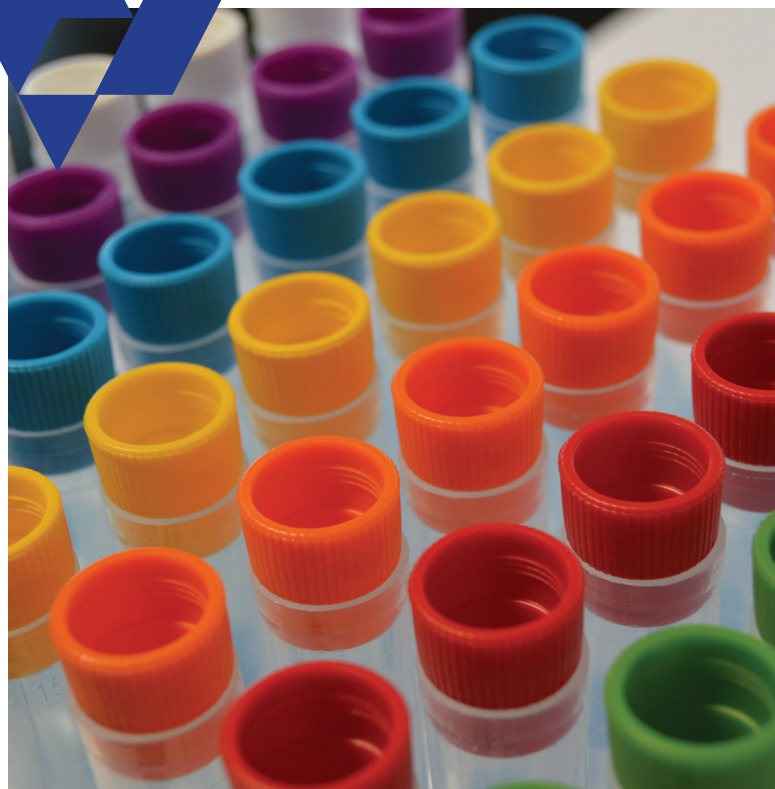

Key Features

- ▶ Stable in the presence of common acids, alkalis and solvents
- ▶ Made with easy-to-read graduation lines that meet Class A tolerances (± 0.10 mL)
- ▶ 8 different colors caps available to help identify samples



Meet the demands of your daily lab routines with AnalytiChem autosampler tubes.

Description	Quantity	Part Number
16 x 100 mm (12 mL) Polypropylene autosampler tubes	1000 tubes	010-515-609
	2 bags of 1000 tubes	010-515-629
	3 bags of 1000 tubes	010-515-639
17 x 100 mm (15 mL) Polypropylene autosampler tubes	1000 tubes	010-515-607
	2 bags of 1000 tubes	010-515-627
	3 bags of 1000 tubes	010-515-637
17 x 100 mm (15 mL) Polypropylene autosampler tubes HandyPAK	45 x 36 tubes	010-515-617

Color caps manufactured from Low-Density Polyethylene (LDPE) with ultra-low metal content.
All caps are made with a double seal flange to guarantee a leak-free seal with the side wall of the tube.

Description	Quantity	Part Number
White	1000 caps	010-515-620
Yellow	1000 caps	010-515-630
Orange	1000 caps	010-515-640
Red	1000 caps	010-515-650
Blue	1000 caps	010-515-660
Clear	1000 caps	010-515-672
Green	1000 caps	010-515-680
Purple	1000 caps	010-515-690

Storage Racks

Description	Quantity	Part Number
Foam Rack for 16x100 autosampler tubes	5	010-515-032
Foam Rack for 17x100 autosampler tubes	5	010-515-038